A EMPRESA
THE COMPANY
LA EMPRESA

OPINIÃO PROFISSIONAL
PROFESSIONAL OPINION
OPINIÓN PROFESIONAL

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THE SYSTHEx® IS A DENTAL IMPLANT COMPANY THAT WAS CREATED TO DEVELOP HIGH QUALITY PRODUCTS AND TO FACILITATE SURGICAL AND PROSTHETIC PROTOCOLS TO CLINICIANS.

THE SYSTHEx® FACTORY POSSESSES THE LATEST TECHNOLOGICAL EQUIPMENTS THAT ARE IN CONFORMITY WITH BOTH NATIONAL AND INTERNATIONAL REGULATIONS. THE SYSTHEx® LABORATORIES OF SURFACE TREATMENT, STERILIZATION, PACKAGING AND LABELING ARE IN AGREEMENT WITH THE HIGHEST STANDARDS FOR PRODUCTION OF IMPLANT PRODUCTS. THE COMPANY HAS THE PRINCIPAL OBJECTIVE TO ATTEND THE MOST EXIGENT CLINICIANS OVER THE WORLD.

SYSTHEx COMPANY IS TESTED BY ANVISA, BY BIOMECHANICAL LABORATORIES (UFSCAR, UNESP) AND MOST IMPORTANTLY, BY CLINICIANS THAT ALREADY ELECTED THE BEST SOLUTION IN IMPLANT DENTISTRY.

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DEAR DOCTOR,

I AM VERY PLEASED TO WELCOME YOU TO THE SYSTHEx IMPLANT SYSTEM. THIS SYSTEM WAS DESIGNED BY A GROUP OF VERY EXPERIENCED CLINICIANS, INCLUDING INTERNATIONAL RESEARCHERS ENGAGED TO PRODUCE EXTREMELY EFFICIENT PRODUCTS, AND TO TRANSFORM IMPLANT DENTISTRY INTO A SIMPLE AND EVEN MORE PROFITABLE ALTERNATIVE. THE COMPANY IS ENGAGED TO CREATE AND TO DEVELOP INNOVATIVE PRODUCTS IN ORDER TO MAKE YOUR PRACTICE SUCCESSFUL. IN ADDITION, OUR RESEARCHES ARE PERFORMED TO BETTER RESPOND TO THE NEEDS OF YOUR PATIENTS. THE SYSTHEx HAVE GREAT SATISFACTION IN BEING ABLE TO HELP YOU.

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RDC 59 - ANVISA - CERTIFICATE OF GOOD MANUFACTURING PRACTICES.

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"THE QUALITY AND PRECISION OF THE SYSTHETX SYSTEM I ALREADY KNEW. FURTHER, THE NEW UNIK® IS AN AMAZING SYSTEM DUE TO ITS STRAIGHTFORWARD CLINICAL PROTOCOL. IN ONLY ONE PACKAGE, THE UNIK PROVIDES WITH ALL THE TOOLS NEEDED TO COMPLETE THE IMPLANT CASE."

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DR. KELLY MISCH

DIPLOMATE, AMERICAN BOARD OF PERIODONTOLGY
FELLOW, INTERNATIONAL CONGRESS OF ORAL IMPLANTOLOGISTS
ADJUNCT PROFESSOR - DEPARTMENT OF PERIODONTICS & ORAL MEDICINE UNIVERSITY OF MICHIGAN - SCHOOL OF DENTISTRY USA

"THE UNIK IMPLANT SYSTEM IS CHARACTERISTIC AND UNIK DUE TO IT'S TREMENDOUS EASE OF USE WITH MINIMAL COMPONENTS. THE HIGH QUALITY AND SIMPLICITY MAKE IT A MASTERFUL PRODUCT FOR BOTH THE IMPLANT SURGEON AND RESTORATIVE DENTIST TO OBTAIN A PRACTICAL, ESTHETIC AND LONG-TERM RESULT."

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"BEYOND THE UNQUESTIONABLE QUALITY FOUND IN ALL THE SYSTHETX SYSTEM, THE UNIK® PRESENTS WITH AN EXCLUSIVE FEATURE: AN EASY PROTOCOL. HENCE, THIS FEATURE GIVES CONFIDENCE AND PREDICTABILITY TO BOTH SURGICAL AND RESTORATIVE PHASES OF IMPLANT THERAPY."

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"THE NEW UNIK SYSTEM BROUGHT AGILITY TO THE LABORATORIAL PHASE OF THE IMPLANT CASE. FOR MANY YEARS, THE SYSTHETX IS THE EXCLUSIVE DISTRIBUTOR OF THE PROSTHETIC COMPONENTS FOR OUR LABORATORIES. WITHOUT ANY QUESTION, THE SYSTHETX IMPRESSES US BY THE HIGH QUALITY AND PRECISION OF ITS PROSTHETIC COMPONENTS."

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DR. DAN ANG

D.D.S., M.S. UNIVERSITY OF MISSOURI KANSAS CITY

“FINALLY A VERSATILE IMPLANT SYSTEM THAT CONSIDERS NOT ONLY THE LONG-TERM BIOLOGICAL EFFECTS OF BONE, BUT ALSO CONSIDERS ESTHETIC RESTORABILITY.”

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DR. DEAN N. AHMAD

DIPLOMATE, AMERICAN BOARD OF PERIODONTOLGY
FELLOW, INTERNATIONAL CONGRESS OF ORAL IMPLANTOLOGY
CLINICAL ASSISTANT PROFESSOR, UNIVERSITY OF CALIFORNIA, SAN FRANCISCO
CHAIRMAN, DENTAL HEALTH COMMITTEE, SACRAMENTO DISTRICT DENTAL SOCIETY
PRIVATE PRACTICE LIMITED TO PERIODONTICS & IMPLANT DENTISTRY (CALIFORNIA, USA).

“AS A PRACTICING PERIODONTIST/IMPLANTOLOGIST IN CALIFORNIA, AND CONSTANTLY BEING FACED WITH HIGH AESTHETIC DEMANDS AND LONG-TERM PREDICTABILITY, THE SYSTHEX IMPLANT SYSTEM HAS PROVED TO BE AN INNOVATIVE TOOL IN MY PRACTICE WHICH MAKES PLACING AND RESTORING DENTAL IMPLANTS WITH EASE.”

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DR. ASVIN VASANTHAN

DIPLOMATE, AMERICAN BOARD OF PERIODONTOLGY
SPECIALTY AND MATER'S OF SCIENCE AT UNIVERSITY OF ALABAMA, BIRMINGHAM.
DIRECTOR OF IMPLANTOLOGY AND CLINICAL ASST. PROFESSOR
UNIVERSITY OF MISSOURI, U.S.A.

“I FEEL THE CONCEPT OF THE UNIK IMPLANT TO BE GREAT. IT IS A WELL DESIGNED IMPLANT WITH A VARIETY OF OPTIONS FOR DIFFERENT CLINICAL SITUATIONS. UNIK IS A GREAT ADDITION TO THE EXISTING IMPLANT MARKET”

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“THE SYSTHEX SYSTEM STUNNED ME BY ITS QUALITY OF ITS SURFACE. I RECOMMEND THIS SYSTEM FOR ALL CASES OF POOR QUALITY BONE”

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DR. TOLGA FİKRET TÖZÜM

ASSOCIATE PROFESSOR, HACETTEPE UNIVERSITY, FACULTY OF DENTISTRY,
DEPARTMENT OF PERIODONTOLGY, TURKEY.
DIPLOMATE, MEFFERT IMPLANT INSTITUTE
FELLOW, INTERNATIONAL CONGRESS OF ORAL IMPLANTOLOGISTS

“I HAVE BEEN INVOLVED IN THE FIELD OF IMPLANTOLOGY FOR MANY YEARS. THE SYSTHEX SYSTEM IS TRULY A PREDICTABLE AND SIMPLE IMPLANT SYSTEM THAT WORKS WELL IN ALL AREAS OF THE MOUTH.”

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
“THE USE OF SYSTHEX ORTHODONTIC IMPLANTS IN THE MODERN ERA DEMONSTRATES THE SIGNIFICANCE OF THIS TOOL IN THE ORTHODONTIC TREATMENT. WITH THE HIGH PREVALENCE OF TOOTH LOSS IN THE MALAR REGION AMONG THE POPULATION, THE SYSTHEX ORTHODONTIC IMPLANT BECOMES EXTREMELY IMPORTANT FOR THE SUCCESS OF THESE CASES.”

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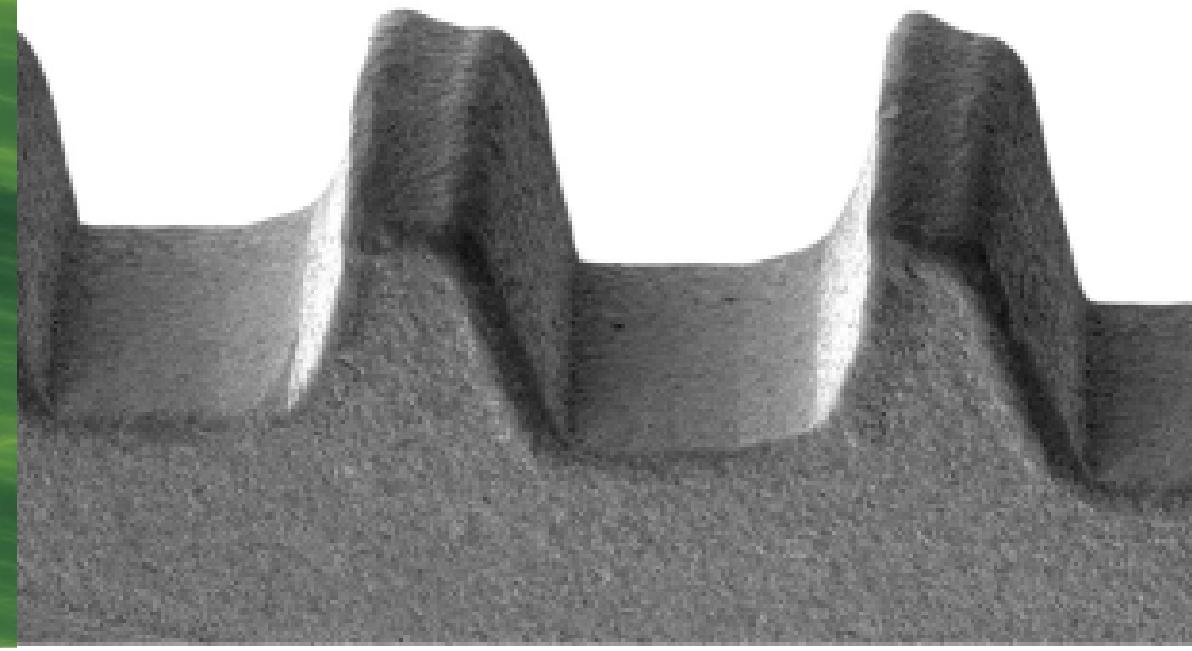The objective of this test was to determine the maximum compression load, based on the Norm ISO 14801: 2003 (Dentistry — Fatigue test for endosseous dental implants).¹

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“THE SYSTHEX PROSTHETIC ABUTMENT DEMONSTRATED TO BE REMARKABLE RESISTANT TO AXIAL FORCES”

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“THE SYSTHEX IMPLANTS DID NOT FRACTURE AFTER MORE THAN FIVE MILLION CYCLES AT 90% OF THE MAXIMUM LOAD.”

O teste foi baseado na Norma ASTM F543:2002 (Standard Specification and Test Methods for Metallic Medical Bone Screws) e NBR ISO 6475-1997.¹

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“THE SYSTHEX IMPLANT PLATFORMS DEMONSTRATED EXCELLENT RESISTANCE WHEN SUBMITTED TO ROTATIONAL FORCES.”

O teste foi baseado na Norma F543:2002 (Standard Specification and Test Methods for Metallic Medical Bone Screws) e NBR ISO 6475-1997.¹

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
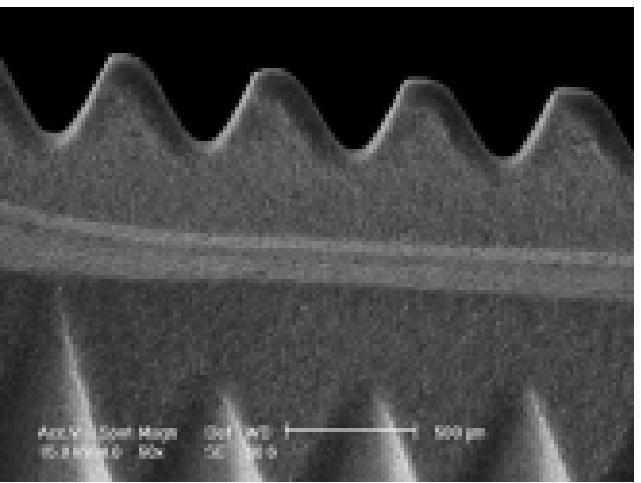
"THE SURFACE OF THE SYSTHEX
IMPLANTS IS TREATED BY SAND
BLASTING + ACID ETCHING,
OFFERING A VERY PREDICTABLE
SURFACE TO ACHIEVE
OSSEointegration."

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SCANNING ELECTRON MICROSCOPE DEMONSTRATING OPTIMAL PORE SIZES. (PERFORMED BY DR. CHARLES
M. COBB, DDS, PhD, PROFESSOR EMERITUS, UNIVERSITY OF MISSOURI-KANSAS CITY, U.S.A.)

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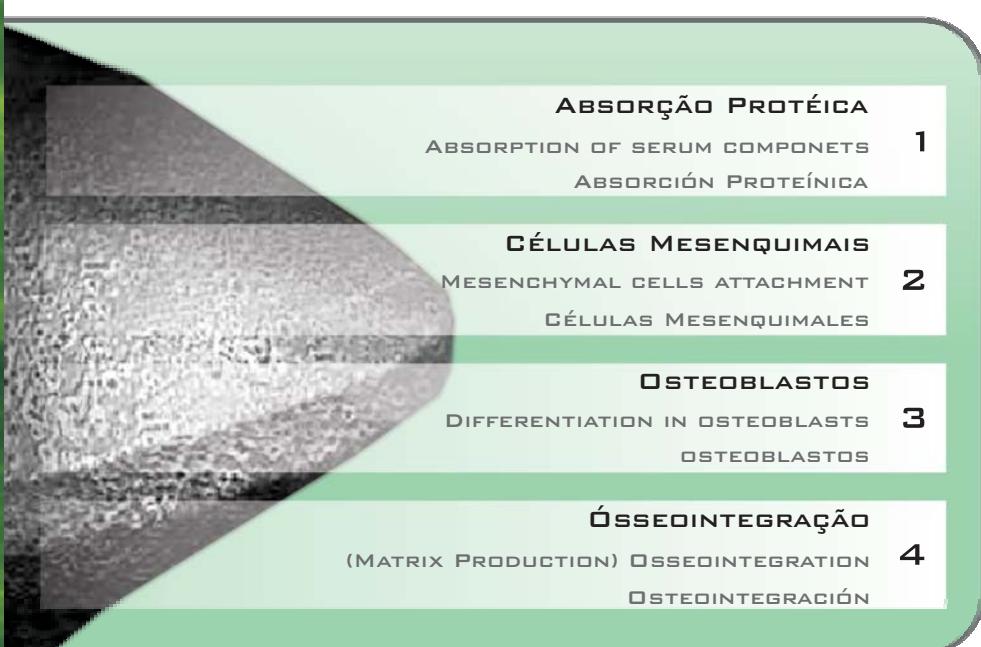
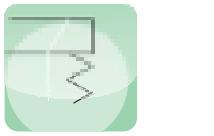
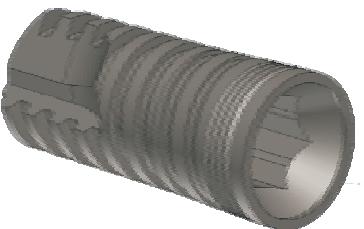
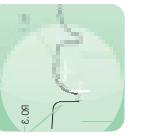
THE UNIK IMPLANT IS AN INNOVATIVE SYSTEM IN WHICH CONTAINS ALL THE NECESSARY COMPONENTS TO PERFORM THE IMPLANT CASE IN ONE SINGLE PACKAGE. THIS ASPECT REPRESENTS NOT ONLY A BETTER COST-BENEFIT TO THE CLINICIAN, BUT ALSO MORE ORGANIZATION TO THE IMPLANT PRACTICE. FURTHER, THE IMPLANT BODY HAS SHALLOW THREADS AT THE CORONAL SIDE TO MINIMIZE THE RISK OF CRESTAL BONE LOSS. ON THE OTHER HAND, THE MID AND APICAL PORTIONS OF THE IMPLANT POSSES DEEP THREADS TO GUARANTEE PRIMARY STABILITY.

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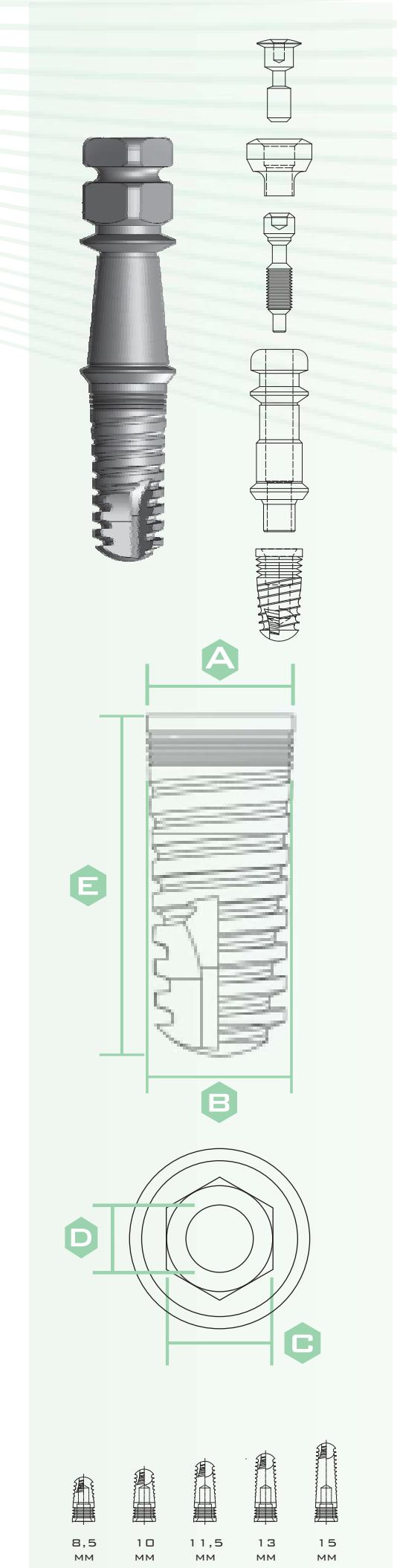COVER SCREW + HEALING
ABUTMENT SCREW IN
ONLY ONE COMPONENT

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CIRCUMFERENTIAL
HEALING ABUTMENT

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PROSTHETIC ABUTMENT
SCREW + TRANSFER SCREW
IN ONLY ONE COMPONENT

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TRANSFER + PROSTHETIC
ABUTMENT IN ONLY ONE COMPONENT

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
SHALLOW THREADS TO MINIMIZE
THE RISK OF CRESTAL BONE LOSS

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DEEP THREADS TO GUARANTEE
PRIMARY STABILITY

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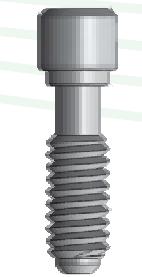BASELINE CLINICAL ASPECT OF TOOTH # 35 (#20) MISSING.

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
CASE DEMONSTRATING ADEQUATE BONE QUANTITY FOR A UNIK IMPLANT PLACEMENT.

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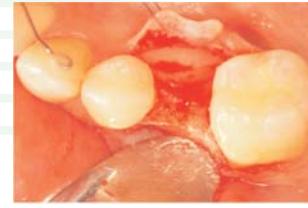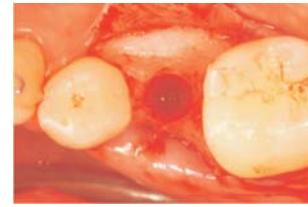
IN THE UNIK PACKAGE, THE TRANSFER + PROSTHETIC ABUTMENT ARE ALREADY INCORPORATED TO THE IMPLANT FOR ITS PLACEMENT.

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UNIK IN ITS FINAL POSITION. AT THIS MOMENT, THE CLINICIAN CAN DECIDE THREE DIFFERENT OPTIONS TO CONTINUE THE CASE:

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2. SINGLE STAGE PROTOCOL WITH HEALING ABUTMENT. THE OPTION CAN BE USED IN CASES WHICH:
A) DO NOT HAVE ENOUGH PRIMARY STABILITY FOR IMMEDIATE LOAD OR FUNCTION (< 30N/cm).
B) THE PROVISIONAL PROSTHESIS AND THE ADJUST OF THE PROSTHETIC ABUTMENT ARE BEING PERFORMED OUT OF THE OFFICE (E.G.: LABORATORY).

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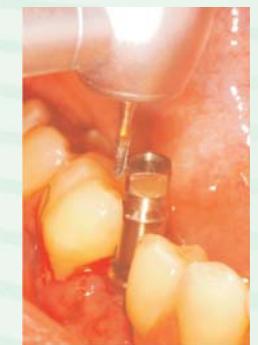
TAKE AN IMPRESSION, PLACE THE TRANSFER (WITH THE MARK AT THE BUCCAL SIDE) AND ANALOG. POOR THE MODEL WITH STONE.

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FOTOGRAFIA DEMONSTRANDO O TRANSFER NO MODELO DE GESSO.

PHOTO DEMONSTRATING THE TRANSFER IN THE DENTAL CAST.

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CUT THE APICAL PORTION OF THE TRANSFER. ADJUST THE PROSTHETIC ABUTMENT ACCORDINGLY.

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
MAKE A PROVISIONAL CROWN AND A METALLIC COPING FOR A FUTURE PICKUP IMPRESSION (PERFORMED WHEN THE IMPLANT IS FINALLY OSSEointegrated).

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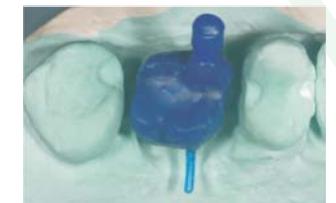PLACE THE PROSTHETIC ABUTMENT WITHOUT APPLYING TORQUE WITH THE WRENCH.

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THIS IMPLANT REPRESENTS THE CLASSIC PRODUCT IN THE FIELD OF IMPLANTOLOGY. "CI" SIGNIFIES INTERNAL TORQUE.

HEXÁGONO EXTERNO

EXTERNAL HEXAGON
STRAIGHT IMPLANT
DOUBLE THREADED
INTERNAL CONNECTION
SURFACE TREATED

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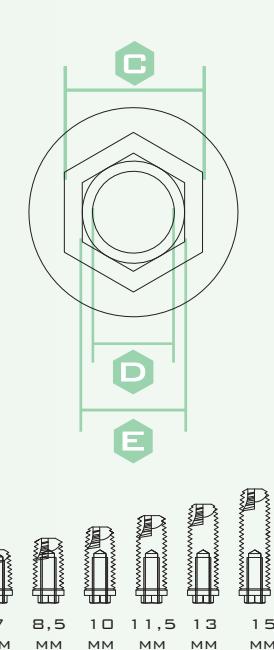
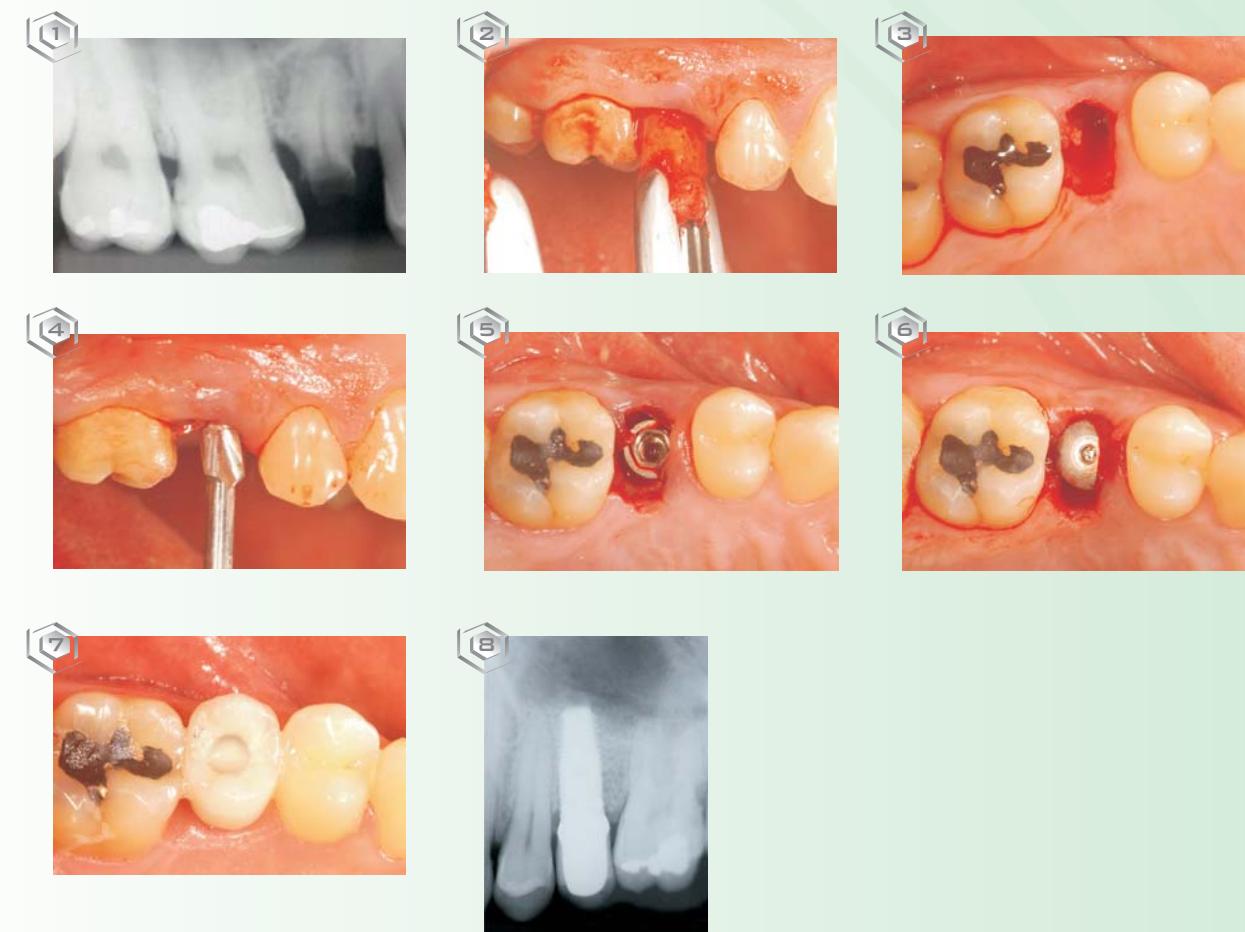
DUPLA ROSCA

CONEXÃO INTERNA

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THE CLASSIC-IN REPRESENTS THE CLASSIC PRODUCT IN THE FIELD OF IMPLANTOLOGY. "IN" SIGNIFIES INTERNAL HEXAGON.

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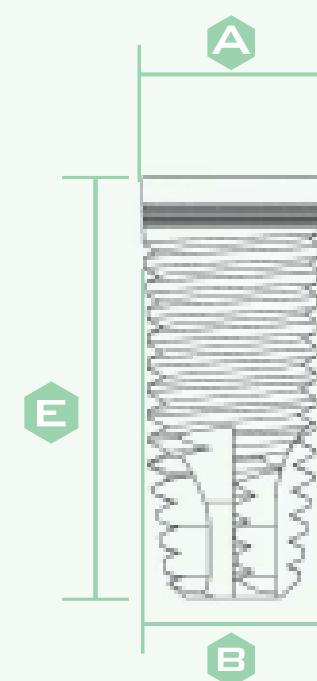
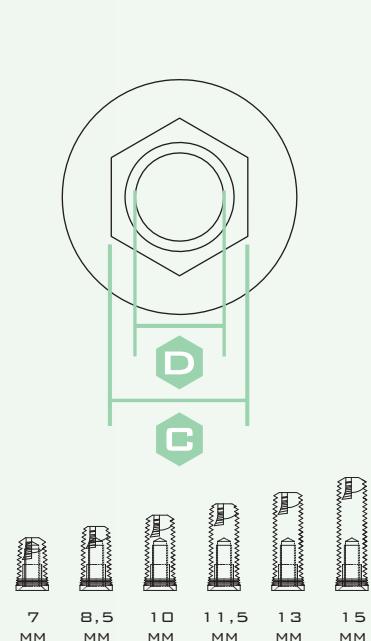
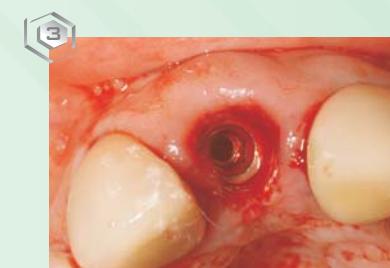
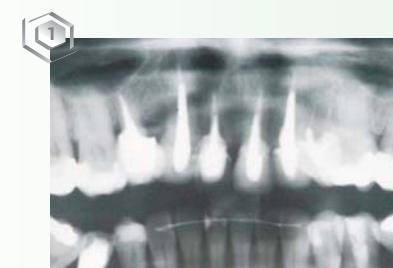
ESTETIC REPRESENTS THE GROUP OF TAPERED IMPLANTS. "CI" SIGNIFIES INTERNAL TORQUE.

EXTERNAL HEXAGON

TAPERED IMPLANT

DOUBLE THREADED

INTERNAL CONNECTION

SURFACE TREATED

DEEP THREADS TO GUARANTEE PRIMARY STABILITY.

FASTER INSTALLATION AND LESS TRAUMA

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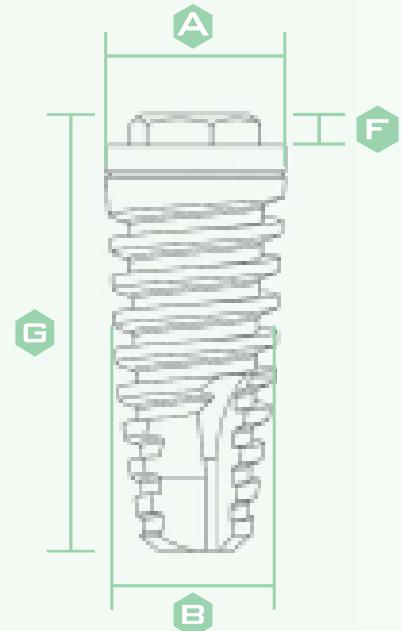
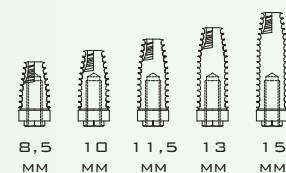
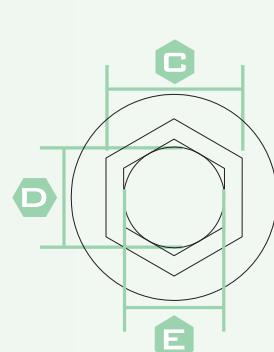
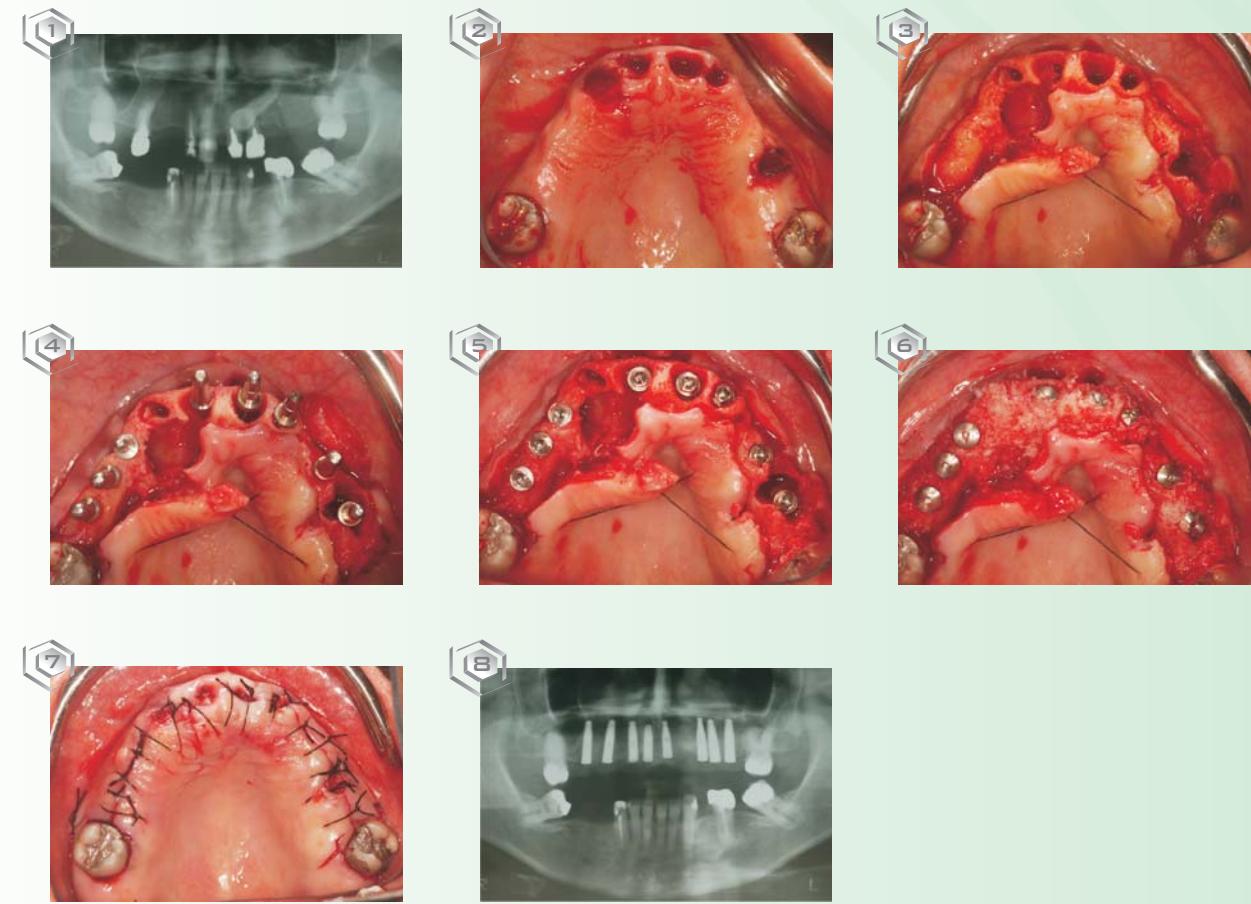
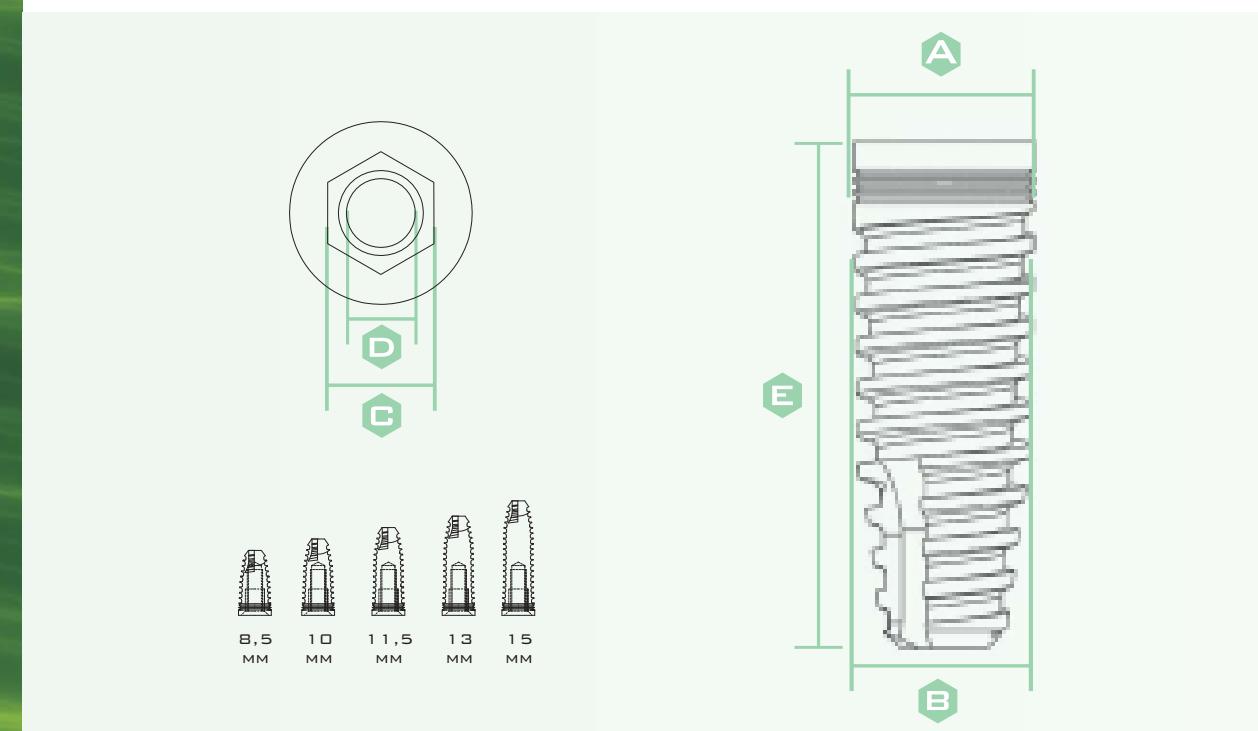
ESTETIC REPRESENTA THE GROUP OF TAPERED IMPLANTS. "IN" SIGNIFIES INTERNAL HEXAGON.

INTERNAL HEXAGON

TAPERED IMPLANT

DOUBLE THREADED

INTERNAL CONNECTION

SURFACE TREATED

DEEP THREADS TO GUARANTEE PRIMARY STABILITY.

FASTER INSTALLATION AND LESS TRAUMA

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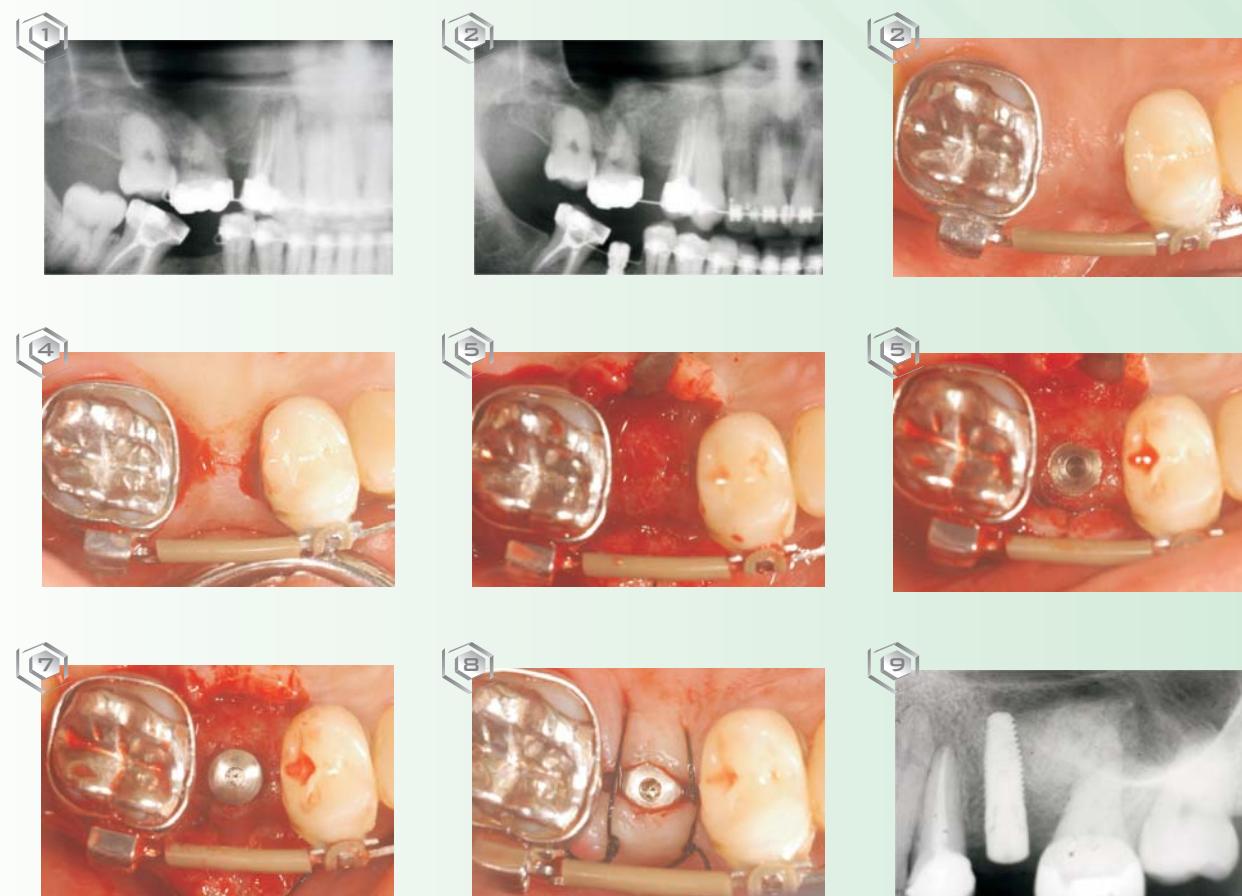
FIT-HEX IS AN IMPLANT WITH 3.3 MM IN DIAMETER. THIS IMPLANT IS DESIGNED TO BE USED IN AREAS OF LIMITED SPACE (SUCH AS LATERAL INCISORS). FIT-HEX IS THE ONLY IMPLANT THAT REQUIRES A MOUNT FOR INSERTION (NOT INCLUDED IN THE SURGICAL KIT).

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1

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PACKAGE PROTECTED WITH THREE LAYERS.

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CHEMICAL INDICATOR BEHIND THE IMPLANT MUST BE IN RED TO INDICATE PROPER GAMA ESTERILIZATION.

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FOR USE, DROP THE LAST PACKAGE IN A STERILE SURFACE.

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WITH STERILE GLOVES, PLACE THE PACKAGE IN A VERTICAL POSITION.

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OPEN THE PACKAGE ONLY WHEN THE SURGICAL PLACEMENT IS CONFIRMED.

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IF YOU DESIRE TO PLACE THE IMPLANT WITH AN ELECTRIC HANDPIECE, USE THE HANDPIECE DRIVER CORRESPONDENT TO THE IMPLANT PLATFORM.

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

IF DESIRED, THE SAME PROTOCOL CAN BE USED WITH MANUAL WRENCH.

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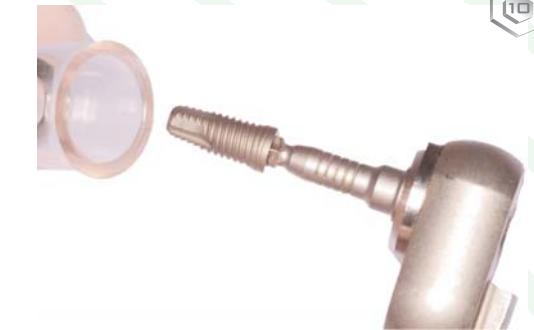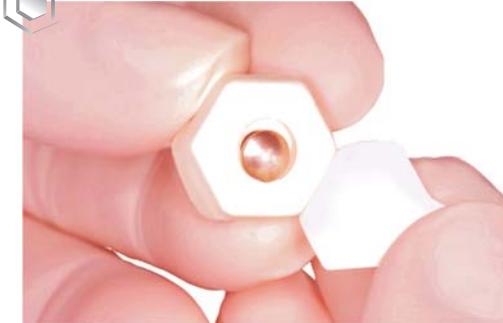THE COVER OF THE PACKAGE CONTAINS THE COVER SCREW PROTECTED WITH A SPECIAL SEAL TO AVOID CONTAMINATION.

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**SINGLE USE
USE BY DATE
BY DATE OF EXPIRATION ONLY
DO NOT USE IF PACKAGE SEAL IS VIOLATED
CONSULT INSTRUCTIONS FOR USE
DO NOT RE-STERILIZE
READ INSTRUCTIONS BEFORE USE**

Description
The package contains one dental implant and one cover screw made of titanium grade 4 (ASTM F67). This material is biocompatible, it allows osseointegration and it is surface treated. The systhex prosthetic system offers both external and internal connections, available in different widths and lengths. The package is protected with three layers (package + blister + primary package).

Indications:
The systhex implants are indicated to replace natural teeth with dental prosthesis. The implants can be used in the upper and lower arches. All dental implants are intended for surgical placement to provide means for pre/prosthetic attachment in single tooth restorations and in partially or fully edentulous spans with multiple single teeth utilizing delayed or immediate loading, or as a terminal or intermediary abutment for fixed or removable bridge work, and to retain overdentures. There are several widths and lengths to be adjusted for each particular case.

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

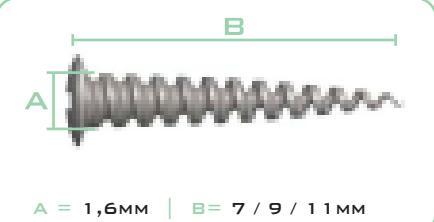
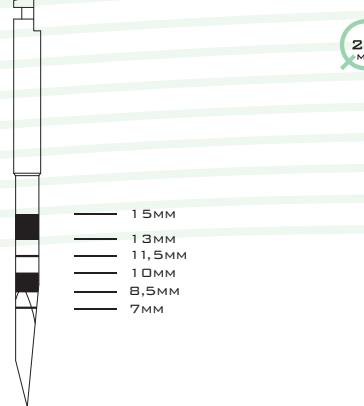
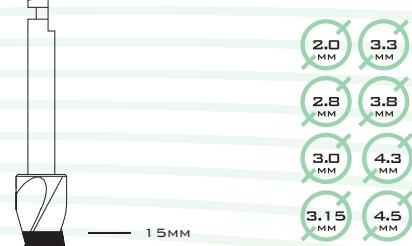
PILOTO
PILOT | PILOTO

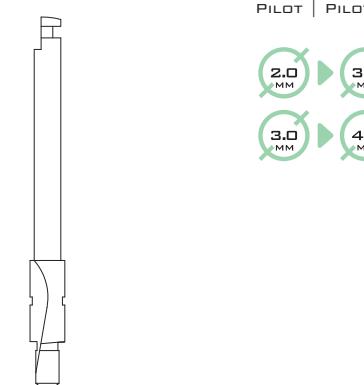PILOT | PILOTO

COUNTERSINKS
COUNTERSINKS | COUNTERSINKS

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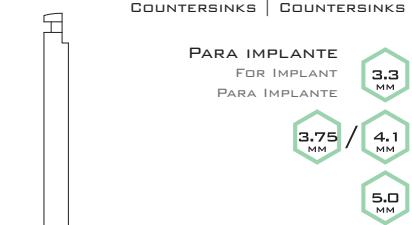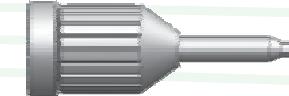
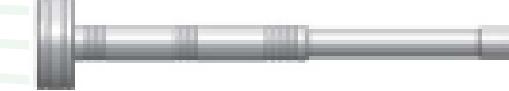
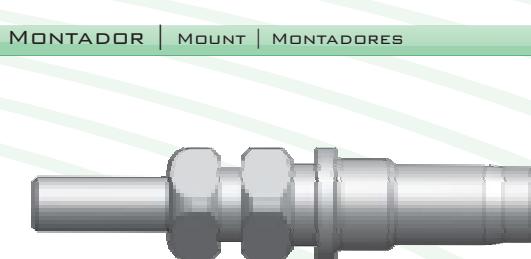
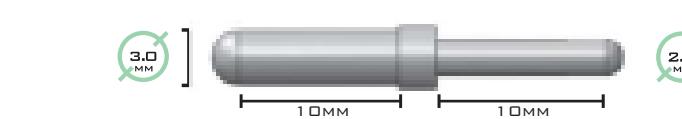
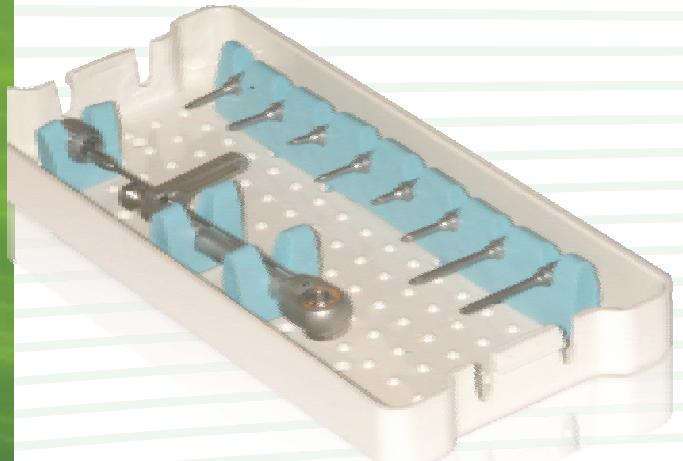
PROSTHETIC KIT | KIT PROTÉSICO

CHAVE TORQUE SEXTAVADA 0.9MM CURTA

0.9 MM HEXAGONAL TORQUE DRIVER SHORT
LLOAVE TORQUE SEXTAVADA 0.9 MM CORTA

CHAVE TORQUE SEXTAVADA 1.2MM CURTA

1.2 MM HEXAGONAL TORQUE DRIVER SHORT
LLOAVE TORQUE SEXTAVADA 1.2 MM CORTA

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
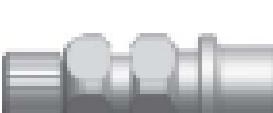
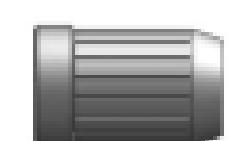
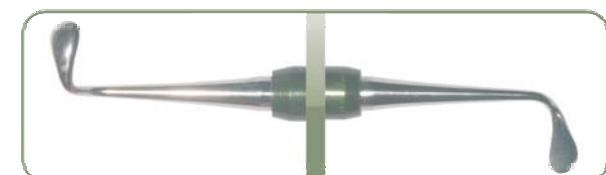
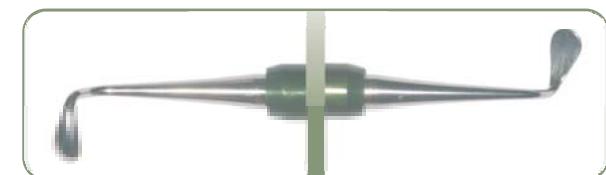
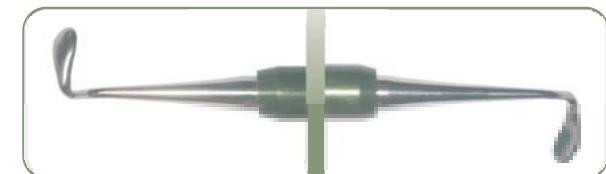
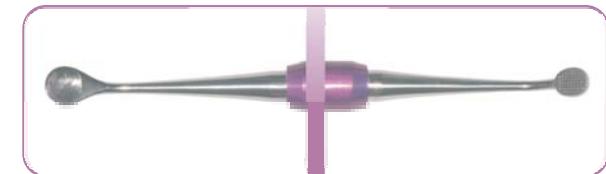
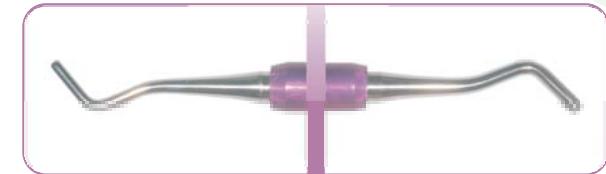
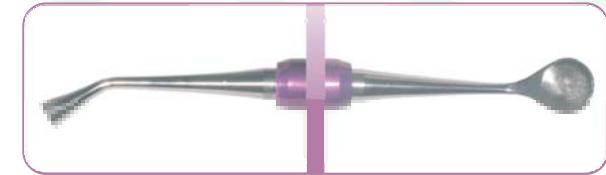
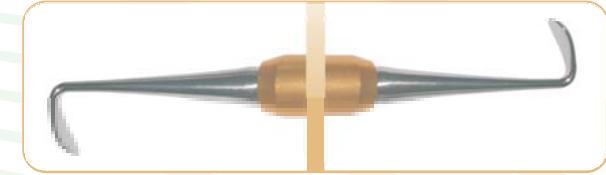
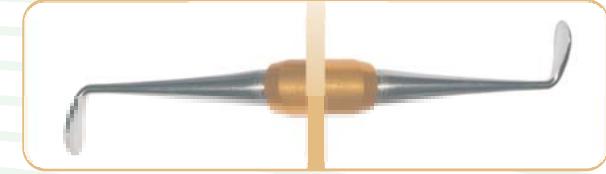
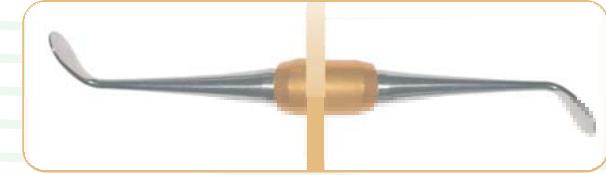
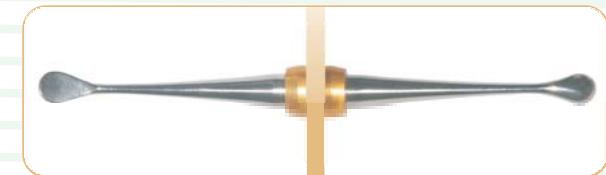
KIT SINUS LIFT
SINUS LIFT KIT | KIT SINUS LIFT

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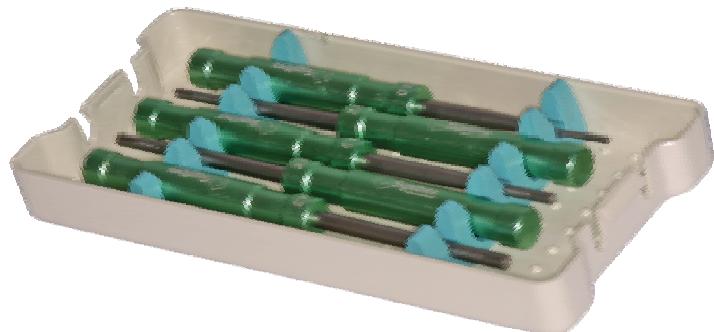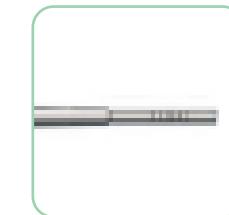
WASHING:
USE NEUTRAL SOAP OR ENZYMAIC DETERGENT
REMOVE DIRT BY HAND OR WITH ULTRASONIC BRUSHES FOR 08 MINUTES.
RINSE WITH DEIONIZED, DISTILLED OR DEMINERALIZED WATER.

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DRYING:
USE SOFT, CLEAN CLOTH OR COMPRESSED AIR

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STERILIZATION:
USE ONLY HEAT STERILIZATION - AUTOCLAVE.
DO NOT STERILIZE IN THE OVEN THE SURGICAL KITS WHOSE MATERIAL IS POLYMER.
NEVER KEEP THE INSTRUMENTAL WET AFTER WASHING.
DRYING THEM PREVENTS OXIDATION.

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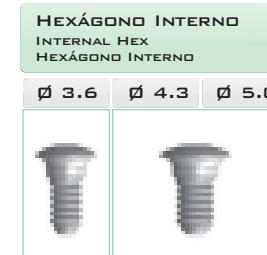
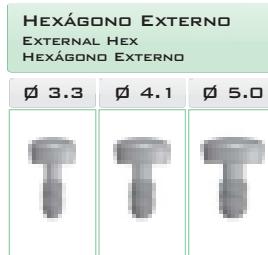
- TITÂNIO.
- TITANIUM
- TITANIO

TRANSFER | TRANSFER | TRANSFER

- DISPONÍVEL EM INOX.
- OPEN OR CLOSED TRAY TECHNIQUE.
- HEXED OR NON-HEXED.
- INCLUDED SCREW.

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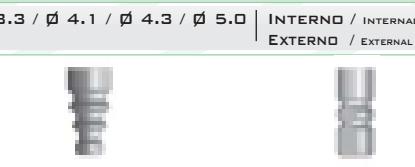
• EXCLUSIVE SYSTHEx SYSTEM:

- ONE COMPONENT ADAPTS TO ALL PLATFORMS.

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
- EXCLUSIVE SYSTHEx SYSTEM:
- ONE COMPONENT ADAPTS TO ALL PLATFORMS
- SIMPLIFIES THE SURGICAL PROCESS

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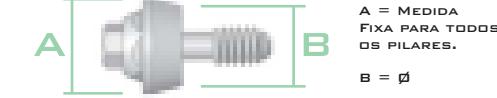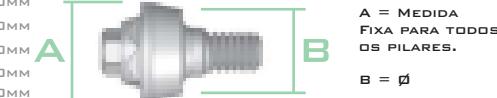
- TITANIUM.
- MINIMAL INTEROCCLUSAL SPACE OF 5.4MM (17°) AND 6.5MM (30°).
- AVAILABLE EXTERNAL HEXAGON
- SCREW INCLUDED

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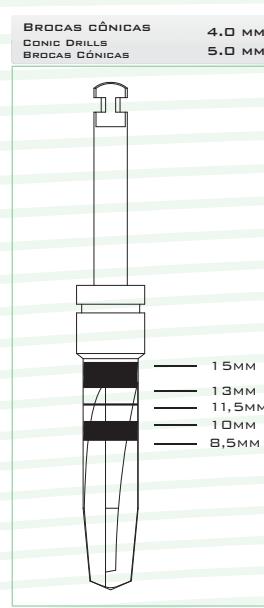
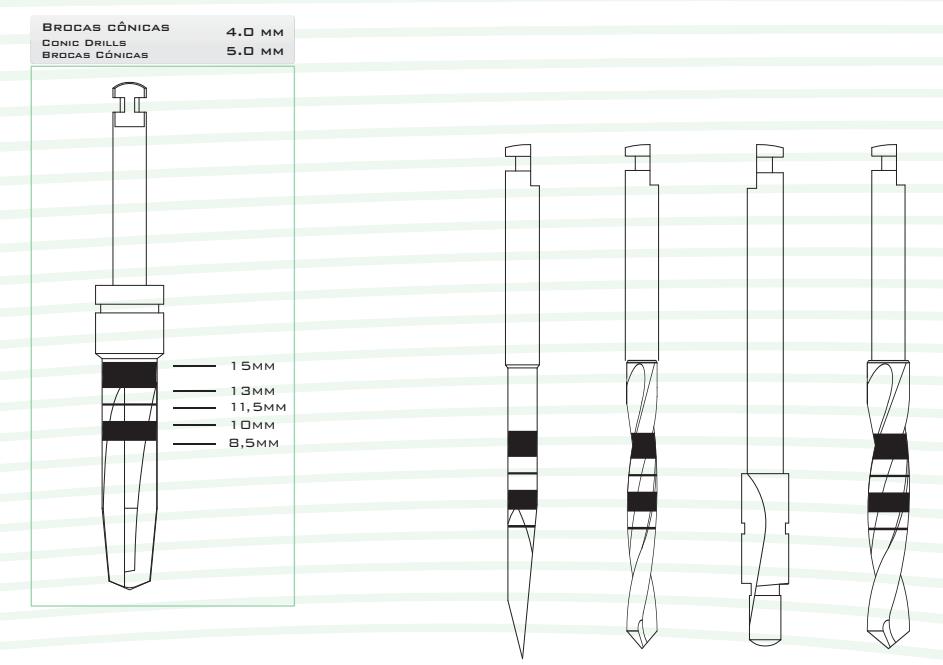
SEQÜÊNCIA DE BROCAS | DRILLING SEQUENCE | SECUENCIA DE BROCAS

3.0 MM	◎	◎	◎	◎	◎	◎	◎	◎	◎	◎	◎	◎	◎	◎
3.15 MM		◎												
COUNTER 3.3 MM			◎											
3.3 MM				◎										
COUNTER 4.1 MM					◎									
3.8 MM						◎								
3/4 MM							◎							
4.3 MM								◎						
4.5 MM									◎					
COUNTER 5 MM										◎				

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THIS DRILL SEQUENCE IS ONLY A SUGGESTION. THE CLINICIAN CAN ALTER THE SEQUENCE ACCORDING TO EACH PARTICULAR CASE.

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814.004	CATRACA CI 5.0 MM CURTA WRENCH DRIVER 5.0 MM SHORT / CARRACA CI 5.0 MM CORTA
814.005	CATRACA CI 5.0 MM MÉDIA WRENCH DRIVER 5.0 MM MEDIUM / CARRACA CI 5.0 MEDIANO
814.006	CATRACA CI 5.0 MM LONGA WRENCH DRIVER 5.0 MM LONG / CARRACA CI 5.0 MM LARGA
814.007	CATRACA IN 3.6 MM CURTA WRENCH DRIVER 3.6 MM SHORT / CARRACA IN 3.6 MM CORTA
814.008	CATRACA IN 3.6 MM MÉDIA WRENCH DRIVER 3.6 MM MEDIUM / CARRACA IN 3.6 MM MEDIANO
814.009	CATRACA IN 3.6 MM LONGA WRENCH DRIVER 3.6 MM LONG / CARRACA IN 3.6 MM LARGA
814.010	CATRACA IN 4.3/5.0 MM CURTA WRENCH DRIVER 4.3/5.0 MM SHORT / CARRACA IN 4.3/5.0 MM CORTA
814.011	CATRACA IN 4.3/5.0 MM MÉDIA WRENCH DRIVER 4.3/5.0 MM MEDIUM / CARRACA IN 4.3/5.0 MM MEDIANO
814.012	CATRACA IN 4.3/5.0 MM LONGA WRENCH DRIVER 4.3/5.0 MM LONG / CARRACA IN 4.3/5.0 MM LARGA
814.013	CONTRA ÂNGULO CI 4.1 MM CONTRA-ANGLE 4.1 MM EXTERNAL HEX / CONTRAÂNGULO CI 4.1 MM
814.014	CONTRA ÂNGULO CI 5.0 MM CONTRA-ANGLE 5.0 MM EXTERNAL HEX / CONTRAÂNGULO CI 5.0 MM
814.015	CONTRA ÂNGULO IN 3.6 MM CONTRA-ANGLE 3.6 MM INTERNAL HEX / CONTRAÂNGULO CI 3.6 MM
814.016	CONTRA ÂNGULO IN 4.3/5.0 MM CONTRA-ANGLE 4.3/5.0 MM INTERNAL HEX / CONTRAÂNGULO CI 4.3/5.0 MM
814.017	CONEXÃO PHILIPS PARA CHAVE DE MÃO PHILLIPS CONNECTOR FOR DRIVER/CONEXION PHILIPS PARA LLAVE MANUAL
814.018	CONEXÃO PHILIPS PARA CONTRA-ÂNGULO PHILLIPS CONNECTOR FOR CONTRA-ANGLE / CONEXION PHILIPS PARA CONTRAÂNGULO
ESTOJOS TRAYS / ESTUCHES	
815.001	KIT CIRÚRGICO (SURGICAL TRAY / KIT QUIRÚRGICO)
815.002	KIT ELEVADORES DE MEMBRANA SINUS LIFT TRAY / KIT LEVANTAMIENTO DE MEMBRANA
815.003	KIT EXPANSORES ÓSSEOS (OSTEOTOME TRAY / KIT EXPANSORES ÓSEOS)
815.004	KIT PROTÉTICO (PROSTHETIC TRAY / KIT DE PRÓTESIS)
815.005	KIT ENXERTO ÓSSEO E ORTODÔNTICO ORTHODONTIC AND GRAFT TRAY / KIT INJERTO ÓSEO Y ORTODÓNTICO
815.006	KIT SINUS LIFT COM 5 CURETAS SINUS LIFT TRAY, 5 CURETTES / KIT SINUS LIFT CON 5 CURETAS
815.007	KIT SINUS LIFT COM 10 CURETAS SINUS LIFT TRAY, 10 CURETTES / KIT SINUS LIFT CON 10 CURETAS
815.100	EXPANSOR ÓSSEO AÇO CIRÚRGICO / ALUMÍNIO 1.8 MM OSTEOTOME STAINLESS STEEL / ALUMINUM / EXPANSOR ÓSSEO ACERO INOXIDABLE / ALUMINIO
815.101	EXPANSOR ÓSSEO AÇO CIRÚRGICO / ALUMÍNIO 2.2 MM OSTEOTOME STAINLESS STEEL / ALUMINUM / EXPANSOR ÓSSEO ACERO INOXIDABLE / ALUMINIO
815.102	EXPANSOR ÓSSEO AÇO CIRÚRGICO / ALUMÍNIO 2.7 MM OSTEOTOME STAINLESS STEEL / ALUMINUM / EXPANSOR ÓSSEO ACERO INOXIDABLE / ALUMINIO
815.103	EXPANSOR ÓSSEO AÇO CIRÚRGICO / ALUMÍNIO 3.0 MM OSTEOTOME STAINLESS STEEL / ALUMINUM / EXPANSOR ÓSSEO ACERO INOXIDABLE / ALUMINIO
815.104	EXPANSOR ÓSSEO AÇO CIRÚRGICO / ALUMÍNIO 3.3 MM OSTEOTOME STAINLESS STEEL / ALUMINUM / EXPANSOR ÓSSEO ACERO INOXIDABLE / ALUMINIO
815.105	ELEVADOR DE MEMBRANA AÇO CIRÚRGICO / ALUMÍNIO 2.0 MM SINUS LIFT STAINLESS STEEL / ALUMINUM / LEVANTADOR DE MEMBRANA ACERO INOXIDABLE / ALUMINIO
815.106	ELEVADOR DE MEMBRANA AÇO CIRÚRGICO / ALUMÍNIO 2.8 MM SINUS LIFT STAINLESS STEEL / ALUMINUM / LEVANTADOR DE MEMBRANA ACERO INOXIDABLE / ALUMINIO
815.107	ELEVADOR DE MEMBRANA AÇO CIRÚRGICO / ALUMÍNIO 3.2 MM SINUS LIFT STAINLESS STEEL / ALUMINUM / LEVANTADOR DE MEMBRANA ACERO INOXIDABLE / ALUMINIO
815.108	ELEVADOR DE MEMBRANA AÇO CIRÚRGICO / ALUMÍNIO 3.8 MM SINUS LIFT STAINLESS STEEL / ALUMINUM / LEVANTADOR DE MEMBRANA ACERO INOXIDABLE / ALUMINIO
815.109	ELEVADOR DE MEMBRANA AÇO CIRÚRGICO / ALUMÍNIO 4.2 MM SINUS LIFT STAINLESS STEEL / ALUMINUM / LEVANTADOR DE MEMBRANA ACERO INOXIDABLE / ALUMINIO
IMPLANTE UNIK UNIK IMPLANT / IMPLANTE UNIK	
816.001	3.75 x 8.5 MM
816.002	3.75 x 10 MM
816.003	3.75 x 11.5 MM
816.004	3.75 x 13 MM
816.005	3.75 x 15 MM
816.006	4.0 x 8.5 MM
816.007	4.0 x 10 MM
816.008	4.0 x 11.5 MM
816.009	4.0 x 13 MM
816.010	4.0 x 15 MM
816.011	5.0 x 8.5 MM
816.012	5.0 x 10 MM
816.013	5.0 x 11.5 MM
816.014	5.0 x 13 MM
816.015	5.0 x 15 MM

IMPLANTE CLASSIC-CI CLASSIC-CI IMPLANT, EXTERNAL HEX / IMPLANTE CLASSIC-CI	
817.001	3.75 x 7 MM
817.002	3.75 x 8.5 MM
817.003	3.75 x 10 MM
817.004	3.75 x 11.5 MM
817.005	3.75 x 13 MM
817.006	3.75 x 15 MM
817.007	4.0 x 7 MM
817.008	4.0 x 8.5 MM
817.009	4.0 x 10 MM
817.010	4.0 x 11.5 MM
817.011	4.0 x 13 MM
817.012	4.0 x 15 MM
817.001C	3.75 x 7 MM
817.002C	3.75 x 8.5 MM
817.003C	3.75 x 10 MM
817.004C	3.75 x 11.5 MM
817.005C	3.75 x 13 MM
817.006C	3.75 x 15 MM
817.007C	4.0 x 7 MM
817.008C	4.0 x 8.5 MM
817.009C	4.0 x 10 MM
817.010C	4.0 x 11.5 MM
817.011C	4.0 x 13 MM
817.012C	4.0 x 15 MM
817.013	4.5 x 7 MM
817.014	4.5 x 8.5 MM
817.015	4.5 x 10 MM
817.016	4.5 x 11.5 MM
817.017	4.5 x 13 MM
817.018	4.5 x 15 MM
817.019	5.0 x 7 MM
817.020	5.0 x 8.5 MM
817.021	5.0 x 10 MM
817.022	5.0 x 11.5 MM
817.023	5.0 x 13 MM
817.024	5.0 x 15 MM
IMPLANTE CLASSIC-IN CLASSIC-IN IMPLANT, INTERNAL HEX / IMPLANTE CLASSIC-IN	
818.001	3.5 x 8.5 MM
818.002	3.5 x 10 MM
818.003	3.5 x 11.5 MM
818.004	3.5 x 13 MM
818.005	3.5 x 15 MM
818.006	3.75 x 7 MM
818.007	3.75 x 8.5 MM
818.008	3.75 x 10 MM
818.009	3.75 x 11.5 MM
818.010	3.75 x 13 MM
818.011	3.75 x 15 MM
818.012	4.0 x 7 MM
818.013	4.0 x 8.5 MM
818.014	4.0 x 10 MM
818.015	4.0 x 11.5 MM
818.016	4.0 x 13 MM
818.017	4.0 x 15 MM
818.018	5.0 x 7 MM
818.019	5.0 x 8.5 MM
818.020	5.0 x 10 MM
818.021	5.0 x 11.5 MM
818.022	5.0 x 13 MM
818.023	5.0 x 15 MM

IMPLANTE ESTETIC-CI ESTETIC-CI IMPLANT, TAPERED, EXTERNAL HEX / IMPLANTE ESTETIC-CI	
819.001	4.0 x 8.5 MM
819.002	4.0 x 10 MM
819.003	4.0 x 11.5 MM
819.004	4.0 x 13 MM
819.005	4.0 x 15 MM
819.006	5.0 x 8.5 MM
819.007	5.0 x 10 MM
819.008	5.0 x 11.5 MM
819.009	5.0 x 13 MM
819.010	5.0 x 15 MM
IMPLANTE ESTETIC-IN ESTETIC-IN IMPLANT, TAPERED, INTERNAL HEX / IMPLANTE ESTETIC-IN	
820.001	4.0 x 8.5 MM
820.002	4.0 x 10 MM
820.003	4.0 x 11.5 MM
820.004	4.0 x 13 MM
820.005	4.0 x 15 MM
820.006	5.0 x 8.5 MM
820.007	5.0 x 10 MM
820.008	5.0 x 11.5 MM
820.009	5.0 x 13 MM
820.010	5.0 x 15 MM
IMPLANTE FIT-HEX FIT-HEX IMPLANT, 3.3 MM / IMPLANTE FIT-HEX	
821.001	3.3 x 8.5 MM
821.002	3.3 x 10 MM
821.003	3.3 x 11.5 MM
821.004	3.3 x 13 MM
821.005	3.3 x 15 MM
IMPLANTE ORTODÔNTICO ORTHODONTIC IMPLANT / IMPLANTE ORTODÔNTICO	
822.001	1.4 x 6 MM - CINTA DE 1 MM (1MM COLLAR / CINTA DE 1MM)
822.002	1.4 x 8 MM - CINTA DE 1 MM (1MM COLLAR / CINTA DE 1MM)
822.003	1.4 x 10 MM - CINTA DE 1 MM (1MM COLLAR / CINTA DE 1MM)
822.004	1.6 x 6 MM - CINTA DE 1 MM (1 MM COLLAR / CINTA DE 1MM)
822.005	1.6 x 8 MM - CINTA DE 1 MM (1 MM COLLAR / CINTA DE 1MM)
822.006	1.6 x 10 MM - CINTA DE 1 MM (1 MM COLLAR / CINTA DE 1MM)
822.007	1.4 x 6 MM - CINTA DE 2 MM (2 MM COLLAR / CINTA DE 2MM)
822.008	1.4 x 8 MM - CINTA DE 2 MM (2 MM COLLAR / CINTA DE 2MM)
822.009	1.4 x 10 MM - CINTA DE 2 MM (2 MM COLLAR / CINTA DE 2MM)
822.010	1.6 x 6 MM - CINTA DE 2 MM (2 MM COLLAR / CINTA DE 2MM)
822.011	1.6 x 8 MM - CINTA DE 2 MM (2 MM COLLAR / CINTA DE 2MM)
822.012	1.6 x 10 MM - CINTA DE 2 MM (2MM COLLAR / CINTA DE 2MM)
KITS KITS / Kits	
823.001	KIT CIRÚRGICO COMPLETO COMPLETE SURGICAL KIT / KIT QUIRÚRGICO COMPLETO
823.002	KIT ELEVADORES DE MEMBRANA
823.003	KIT SINUS LIFT / KIT LEVANTAMIENTO DE MEMBRANA KIT EXPANSORES ÓSEOS OSTEOTOME KIT / KIT EXPANSORES ÓSEOS
823.004	KIT PROTÉTICO COMPLETO COMPLETE PROSTHETIC KIT / KIT PRÓTESIS COMPLETO
823.005	KIT SINUS LIFT COM 5 CURETAS SURGICAL KIT 5 CURETTES / KIT SINUS LIFT CON 5 CURETAS
823.006	KIT SINUS LIFT COM 10 CURETAS EM ALUMÍNIO SURGICAL KIT 10 ALUMINUM CURETTES / KIT SINUS LIFT CON 10 CURETAS EN ALUMINIO
823.007	KIT SINUS LIFT COM 10 CURETAS EM AÇO INOX SURGICAL KIT 10 INOX CURETTES / KIT SINUS LIFT CON 10 CURETAS EN ACERO INOX
823.008	KIT ENXERTO ÓSSEO E ORTODÔNTICO BONE GRAFT AND ORTHODONTIC ANCHORAGE TRAY KIT INJERTO ENXERTO ÓSEO Y ORTODÓNTICO
823.009	KIT ENXERTO ÓSSEO BONE GRAFT ANCHORAGE TRAY / KIT INJERTO ENXERTO ÓSEO
823.010	KIT ORTODÔNTICO ORTHODONTIC ANCHORAGE TRAY / KIT ORTODÔNTICO
MACHO ROSCA TAP DRILL / MACHO ROSCA	
824.001	3.75 MM
824.002	4.0 MM
824.003	5.0 MM

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